

User Guide

This guidance should be read in conjunction with the IFU supplied with the specific Breathing Circuit, the Flow Driver and Heater Humidifier

Step 1

Size the patient for prong or mask.




Step 2

Fit the indicated colour coded nCPAP prong or mask to the nCPAP generator.



Step 3

Fit nCPAP generator pressure and flow lines to NeoFlow® circuit using dedicated connector provided.

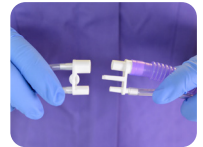

Step 4

To fit support cradle, approximate distance from nares to brow.



Step 5

Attach support cradle to nCPAP generator. Ensure pressure and gas supply lines are fitted in slot provided and exhaust line is secured above.



Step 6

Apply nCPAP prong or mask and secure support cradle to bonnet.



Step 7

Apply the grey securement lines to Velcro® tabs on the bonnet. Ensure prong or mask is securely fitted in nares but **NOT** overtight.

